

LEGEND MAX™

ELISA Kit with Pre-coated Plate



Human IGFALS

Cat. No. 445907

ELISA Kit for Accurate Quantitation of Human IGFALS from Cell Culture Supernatant, Serum, Plasma and Other Biological Fluids

BioLegend, Inc. biolegend.com

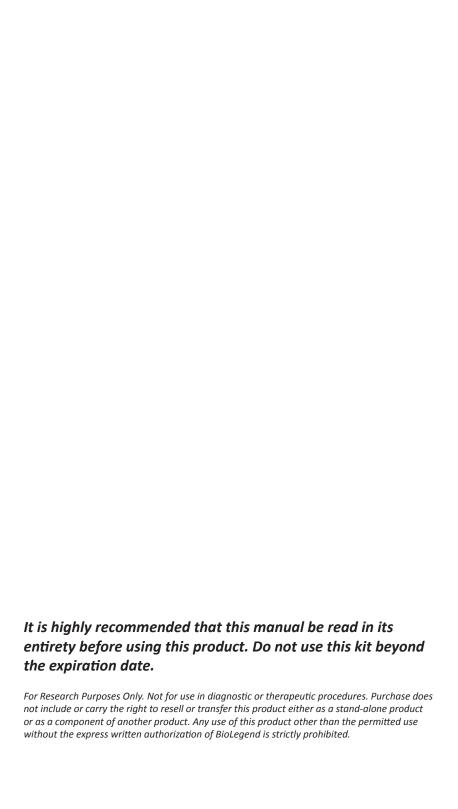




Table of Contents	Page
Introduction	2
Materials Provided	. 2
Materials to be Provided by the End-User	. 3
Storage Information	3
Health Hazard Warnings	. 3
Specimen Collection and Handling	· 4
Reagent and Sample Preparation	4
Assay Procedure	5
Assay Procedure Summary	. 7
Calculation of Results	. 8
Typical Data	8
Performance Characteristics	. 9
Specificity	. 9
Sensitivity	. 9
Recovery	9
Linearity	9
Intra-Assay Precision	10
Inter-Assay Precision	10
Biological Samples	10
Troubleshooting Guide	11
ELISA Plate Template	. 12

Introduction:

Insulin-like growth factor (IGF) is an important hormone for childhood growth and in a variety of physiological functions in adults. Free IGF has a short half-life in serum of about 10 minutes. This half-life is extended by forming a binary complex with IGFBP1-6 to 30 – 90 minutes. IGFALS binds the IGF/IGFBP3 or IGBP5 complex to generate a ternary complex with a serum half-life of more than 12 hours. About 80% of circulating IGF is in this ternary complex.

Insulin-like growth factor binding protein, acid labile subunit (IGFALS), is produced by the liver and is present at high concentrations in postnatal serum. IGFALS is encoded by the IGFALS gene and stimulated by growth hormone. Defects in the IGFALS gene result in lower circulating levels of IGF-1 and IGF-2, causing delayed puberty, increased risk of metastatic disease, and is associated with insulin resistance. Different conditions such as fasting, catabolic disease, or cirrhosis have been shown to reduce serum IGFALS levels in humans. The IGFALS protein is composed of 578 amino acids with 20 leucine-rich repeats.

BioLegend's LEGEND MAX™ Human IGFALS ELISA kit is a Sandwich Enzyme-Linked Immunosorbent Assay (ELISA) with a 96-well strip plate that is precoated with a monoclonal mouse anti-human IGFALS capture antibody. The detection antibody is a biotinylated monoclonal mouse anti-human IGFALS antibody. This kit is specifically designed for the accurate quantitation of human IGFALS from cell culture supernatant, serum, plasma and other biological fluids. Free or ternary complex IGFALS can be detected. It is analytically validated with ready-to-use reagents.

Materials Provided:

Description	Quantity	Volume	Part #
Anti-Human IGFALS Pre-coated 96- well Strip Microplate	1 plate		75000
Human IGFALS Detection Antibody	1 bottle	12 mL	75001
Human IGFALS Standard	1 vial	lyophilized	75003
Avidin-HRP	1 bottle	12 mL	77897
Assay Buffer B	1 bottle	25 mL	79128
Wash Buffer (20X)	1 bottle	50 mL	78233
Substrate Solution F	1 bottle	12 mL	76335
Stop Solution	1 bottle	12 mL	79133
Plate Sealers	4 sheets		78101

Materials to be Provided by the End-User:

- Microplate reader able to measure absorbance at 450 nm
- Adjustable pipettes to measure volumes ranging from 1 μL to 1,000 μL
- · Deionized water
- Wash bottle or automated microplate washer
- Log-Log graph paper or software for data analysis
- Polypropylene tubes to prepare standard dilutions
- Timer
- Plate Shaker
- Polypropylene vials

Storage Information:

Store unopened kit components between 2°C and 8°C. Do not use this kit beyond its expiration date.

O	Opened or Reconstituted Components				
Microplate wells	If not all microplate strips are used, remove the excess strips by pressing up from underneath each strip. Place excess strips back in the foil pouch with the included desiccant pack and reseal. Store between 2°C and 8°C for up to one month.				
Standard	The remaining reconstituted standard stock solution can be aliquoted into polypropylene vials and stored at -70°C for up to one month. Avoid repeated freeze-thaw cycles.				
Detection Antibody					
Avidin-HRP					
Assay Buffer B	Store opened reagents between 2°C and 8°C and use within one month.				
Wash Buffer (20X)					
Substrate Solution F					
Stop Solution					

Health Hazard Warnings:

- Reagents that contain preservatives may be harmful if ingested, inhaled or absorbed through the skin. Refer to the MSDS online at BioLegend's website for details (www.biolegend.com/msds).
- 2. Substrate Solution F is harmful if inhaled or ingested. Avoid skin, eye and clothing contact.

- 3. To reduce the likelihood of blood-borne transmission of infectious agents, handle all serum, plasma and other biological fluids in accordance with NCCLS regulations.
- 4. Stop Solution contains strong acid. *Wear eye, hand, and face protection.*
- 5. Before disposing of the plate, rinse it with an excess amount of tap water.

Specimen Collection and Handling:

Specimens should be clear and non-hemolyzed. If possible, unknown samples should be run at a number of dilutions to determine the optimal dilution factor that will ensure accurate quantitation.

<u>Cell Culture Supernatant</u>: If necessary, centrifuge all samples to remove debris prior to analysis. It is recommended that samples be stored at < -70°C. Avoid repeated freeze-thaw cycles.

<u>Serum</u>: Use a serum separator tube and allow clotting for at least 30 minutes, then centrifuge for 10 minutes at 1,000 x g. Remove serum layer and assay immediately or store serum samples at < -70°C. Avoid repeated freeze-thaw cycles.

<u>Plasma:</u> Collect blood samples in citrate, heparin or EDTA containing tubes. Centrifuge for 10 minutes at 1,000 x g within 30 minutes of collection. Assay immediately or store plasma samples at < -70°C. Avoid repeated freeze-thaw cycles.

Reagent and Sample Preparation:

Note: All reagents should be diluted immediately prior to use.


- Dilute the 20X Wash Buffer to 1X with deionized water. For example, make 1 liter of 1X Wash Buffer by adding 50 mL of 20X Wash Buffer to 950 mL of deionized water. If crystals have formed in the 20X Wash Buffer, bring to room temperature and vortex until dissolved.
- Reconstitute the lyophilized Human IGFALS Standard by adding the volume of Assay Buffer B to make the 480 ng/mL standard stock solution (Refer to LEGEND MAX Kit Lot-Specific Certificate of Analysis/LEGEND MAX Kit Protocol). Allow the reconstituted standard to sit at room temperature for 15-20 minutes, then briefly vortex to mix completely.
- 3. In general, serum and plasma samples require 80-fold dilution and tissue culture supernatant do not require dilution. If dilution is required, use Assay Buffer B. Lower or higher dilutions can be used as needed, but must be quantified to fit the assay range by the end user.

⁴ Tel: 858-768-5800

Assay Procedure:

Note: Do not mix reagents from different kits or lots. Reagents and/or antibodies from different manufacturers should not be used with this kit.

- 1. Bring all reagents to room temperature prior to use. It is strongly recommended that all standards and samples be run in duplicate or triplicate. A standard curve is required for each assay.
- 2. If not all microplate strips will be used, remove the excess strips by pressing up from underneath each strip. Place excess strips back in the foil pouch with the included desiccant pack and reseal.
- 3. Transfer 400 μL of the 480 ng/mL top standard to a separate test tube. Perform five two-fold serial dilutions of the 480 ng/mL top standard using Assay Buffer B as the diluent. Thus, the Human IGFALS standard concentrations in the tubes are 480 ng/mL, 240 ng/mL, 120 ng/mL, 60 ng/mL, 30 ng/mL, and 15 ng/mL, respectively. Assay Buffer B serves as the zero standard (0 ng/mL).

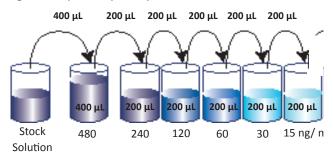


- 4. Wash the plate 4 times with at least 300 μL of 1X Wash Buffer per well and blot any residual buffer by firmly tapping the plate upside down on absorbent paper. All subsequent washes should be performed similarly.
- 5. Add 50 μ L of Assay Buffer B to each well that will contain either standard dilutions or samples. Then add 50 μ L of standard dilutions or samples to the appropriate wells.
- 6. Seal the plate with a Plate Sealer included in the kit and incubate the plate for 2 hours at room temperature with shaking.
- 7. Discard the contents of the plate into a sink, then wash the plate 4 times with 1X Wash Buffer as in step 4.
- 8. Add 100 μ L of Human IGFALS Detection Antibody solution to each well, seal the plate and incubate at room temperature for 1 hour while shaking.
- 9. Discard the contents of the plate into a sink, then wash the plate 4 times


- with 1X Wash Buffer as in step 4.
- 10. Add 100 μ L of Avidin-HRP solution to each well, seal the plate and incubate at room temperature for 30 minutes while shaking.
- 11. Discard the contents of the plate into a sink, then wash the plate 5 times with 1X Wash Buffer as in step 4. For this final wash, soak wells in 1X Wash Buffer for 30 seconds to 1 minute for each wash. This will help minimize background.
- 12. Add 100 μ L of Substrate Solution F to each well and incubate for 15 minutes in the dark. Wells containing human IGFALS should turn blue in color with an intensity proportional to its concentration. It is not necessary to seal the plate during this step.
- 13. Stop the reaction by adding 100 μ L of Stop Solution to each well. The solution color should change from blue to yellow.
- 14. Read absorbance at 450 nm within 30 minutes. If the reader is capable of reading at 570 nm, the absorbance at 570 nm can be subtracted from the absorbance at 450 nm.

Assay Procedure Summary


1. Wash 4 times. Add 50 μL Assay Buffer B to standard wells and sample wells

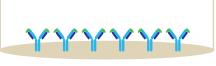


2. Add 50 μL of standard or sample, incubate 2 hr, RT, shaking

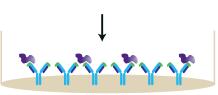


Wash 4 times
 Add 100 µL Detection Antibody solution
 Incubate 1 hr, RT, shaking

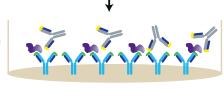

 Wash 4 times Add 100 μL Avidin-HRP solution Incubate 30 min, RT, shaking

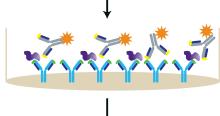

5. Wash 5 times Add 100 μL Substrate Solution F Incubate 15 min, RT, in the dark

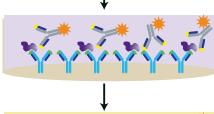

6. Add 100 μL Stop Solution



7. Read absorbance at 450 nm and 570 nm


Calculation of Results:

The data can be best calculated with computer-based curve-fitting software using a 5- or 4-parameter logistics curve-fitting algorithm. If an appropriate software is not available, use log-log graph paper to determine sample concentrations. Determine the mean absorbance for each set of duplicate or triplicate standards, controls, and samples. Plot the standard curve on log-log graph paper with cytokine concentration on the X-axis and absorbance on the Y-axis. Draw a best fit line through the standard points. To determine the unknown cytokine concentrations, find the mean absorbance value of the unknown concentration on the Y-axis and draw a horizontal line to the standard curve. At the point of intersection, draw a vertical line to the X-axis and read the cytokine concentration.

If samples were diluted, multiply the concentration by the appropriate dilution factor. If a test sample's absorbance value falls outside the linear portion of the standard curve, the test sample needs to be re-analyzed at a higher (or lower) dilution as appropriate.

Typical Data:

This standard curve was generated at BioLegend for demonstration purposes only. A standard curve must be run with each assay.

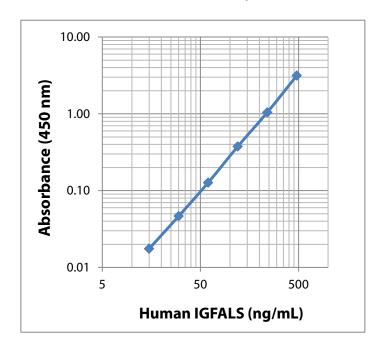

Performance Characteristics:

<u>Specificity:</u> This kit recognizes natural and recombinant human IGFALS. No cross reactivity was observed when this kit was used to analyze the following recombinant proteins at 50 ng/mL.

Human	IGFBP4, CD14, IGFBP3, IGFBP6, IGFBP7, IGF-1, IGF-2, MMP1,
	EGF, EPO, FASL, FGF-basic, GMCSF, IFNg, IL-3, IL-4, Lta (TNFb),
	MCSF, NGAL, IL-6, IL-8, IL-21, IL-22, IL-23, IL-27, IL-33, IL-34,
	Paxillin, PLGF-1, SCF, TGFa, TIMP-1, TNFa, Trance (Rankl)
Mouse	CCL5, CD74, IGFBP1, IGFBP2, IGFBP4, IGFPB5, IGFBP6, IGF-1,
	IGF-2

No reactivity was detected with mouse, rat, bovine, or rabbit serum. Slight reactivity detected with goat serum. Strong reactivity detected with porcine serum.

<u>Interference:</u> Ternary complexes with IGFALS/IGFBP3/IGF-1 and IGFALS/IGFBP3/IGF-2 were assayed. No interference from complex formation was detected.

<u>Sensitivity:</u> The average minimum detectable concentration of human IGFALS is 4.39 ± 0.57 ng/mL (n=10).

<u>Recovery:</u> Recombinant human IGFALS (240, 120 and 60 ng/mL) was spiked into human samples, then analyzed with the LEGEND MAX™ Human IGFALS ELISA kit.

Sample Type	N	% Recovery
Serum	5	106.9
Citrate Plasma	5	107.0
EDTA Plasma	5	99.5
Heparin Plasma	5	105.0

<u>Linearity:</u> Human samples containing high levels of IGFALS were diluted to produce samples within the dynamic range of the kit, and then assayed to determine the dilutional linearity.

Sample Type	N	% Linearity
Serum	5	108.2
Citrate Plasma	5	104.9
EDTA Plasma	5	114.9
Heparin Plasma	5	115.5

<u>Intra-Assay Precision:</u> Two human samples containing different IGFALS concentrations were tested on one plate with 16 replicates.

Concentration	Sample 1	Sample 2
Number of Replicates	16	16
Mean Concentration (pg/mL)	215.2	42.2
Standard Deviation	18.3	3.6
%CV	8.5	8.4

<u>Inter-Assay Precision:</u> Two human samples containing different IGFALS concentrations were tested in four independent assays

Concentration	Sample 1	Sample 2
Number of Assays	4	4
Mean Concentration (pg/mL)	225.3	55.3
Standard Deviation	11.5	5.3
%CV	5.1	7.8

Normal Biological Samples:

Normal Ranges: Human serum and plasma samples (n = 32) were assayed for natural human IGFALS. Serum and plasma range detected was $7.3-13.9\ \mu g/mL.$ No significant difference was detected between serum, citrate plasma, EDTA plasma, and Heparin plasma in paired samples.

Troubleshooting Guide:

Problem	Probable Cause	Solution		
High Background	Background wells were contaminated	Avoid cross-well contamination by using the provided plate sealers. Use multichannel pipettes and change tips between pipetting samples and reagents.		
	Insufficient washes	Increase number of washes. Increase soaking time between washes prior to addition of substrate solution.		
	TMB Substrate Solution was contaminated	TMB Substrate Solution should be clear and colorless prior to addition to wells. Use a clean container prior to pipetting substrate solution into wells.		
No or poor signal	Detection Antibody, Avidin-HRP or Substrate solution were NOT added			
	Wrong reagent or reagents were added in wrong sequential order	Rerun the assay and follow the protocol.		
	Insufficient plate agitation	The plate should be agitated during all incubation steps using a plate shaker at a speed where solutions in wells are within constant motion without splashing.		
	The wash buffer contains Sodium Azide (NaN3)	Avoid Sodium Azide contamination in the wash buffer as it inhibits HRP activity.		
	Incubations were done at an inappropriate temperature, timing or without agitation	Rerun the assay and follow the protocol.		
Low or poor standard curve	The standard was incorrectly reconstituted or diluted	Adjust the calculations and follow the protocol.		
signal	Standard was inappropriately stored	Store the reconstituted standard stock solution in polypropylene vials at -70°C. Avoid repeated freeze-thaw cycles.		
	Reagents added to wells with incorrect concentrations	Check for pipetting errors and the correct reagent volume.		

Problem	Probable Cause	Solution		
Signal is high, standard curves have saturated	Standard reconstituted with less volume than required	Reconstitute new lyophilized standard with the correct volume of solution recommended in the protocol.		
signal	Standards/samples, detection antibody, Avidin-HRP or substrate solution were incubated for too long	Rerun the assay and follow the protocol.		
Sample readings	Samples contain no or below detectable levels of the analyte	If samples are below detectable levels, it may be possible to use a larger sample volume. Contact technical support for appropriate protocol modifications.		
are out of range	Samples contain analyte concentrations greater than highest standard point	Samples may require dilution and analysis.		
	Multichannel pipette errors	Confirm that pipette calibrations are accurate.		
High variation in	Plate washing was not adequate or uniform	Ensure pipette tips are tightly secured.		
samples and/or standards	Non-homogenous samples	Ensure uniformity in all wash steps. Thoroughly mix samples before assaying.		
	Samples may have high particulate matter	Remove particulate matter by centrifugation.		
	Cross-well contamination	Do not reuse plate sealers.		
		Always change tips for reagent additions. Ensure that pipette tips do not touch the reagents on the plate.		

	12								
	11								
	10								
	6								
a	8								
ELISA Plate Template	7								
late Te	9								
LISA P	5								
	4								
	3								
	2								
	1								
		A	В	C	D	田	F	Ð	Н



LEGEND MAX™ Kits are manufactured by **BioLegend Inc.**

9727 Pacific Heights Blvd.

San Diego, CA 92121 Tel: 1.858.768.5800

Tel: US & Canada Toll-Free: 1.877.Bio-Legend (1.877.246.5343)

Fax: 1.877.455.9587

Email: info@biolegend.com

biolegend.com

For a complete list of world-wide BioLegend offices and distributors, please visit our website at: biolegend.com