

FITC anti-human CD11c Antibody

Catalog# / Size	301603 / 25 tests 301604 / 100 tests
Clone	3.9
Regulatory Status	RUO
Workshop	III NL707
Other Names	Integrin α X subunit, CR4, p150, ITGAX
Isotype	Mouse IgG1, κ
Description	CD11c is a 145-150 kD type I transmembrane glycoprotein also known as integrin α X and CR4. CD11c non-covalently associates with integrin β 2 (CD18) and is expressed on monocytes/macrophages, dendritic cells, granulocytes, NK cells, and subsets of T and B cells. CD11c has been reported to play a role in adhesion and CTL killing through its interactions with fibrinogen, CD54, and iC3b.

Product Details

Reactivity	Human, African Green, Baboon, Chimpanzee, Cynomolgus, Rhesus, Squirrel Monkey
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography, and conjugated with FITC under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 μ l per million cells in 100 μ l staining volume or 5 μ l per 100 μ l of whole blood.
Excitation Laser	Blue Laser (488 nm)
Application Notes	Clone 3.9 preferentially binds the activated form of CD11c, is specific for the I domain of CD11c, and is able to partially block the binding of CD11c and ICAM-4. 3.9 binding is divalent cation dependent ¹² . While analyzing blood, it is best to use heparin as the anti-coagulant and not EDTA. Since the ability of clone 3.9 to bind to its target is divalent cation dependent, the usage of EDTA as an anti-coagulant may be detrimental to staining due to its chelating properties. Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections ⁴ , and functional assays ^{5,6} . The LEAF™ purified antibody (Endotoxin <0.1 EU/ μ g, Azide-Free, 0.2 μ m filtered) is recommended for functional assays (Cat. No. 301616). For highly sensitive assays, we recommend Ultra-LEAF™ purified antibody (Cat. No. 301632) with a lower endotoxin limit than standard LEAF™ purified antibodies (Endotoxin <0.01 EU/ μ g).

Application References

- Schlossman S, *et al.* Eds. 1995. Leucocyte Typing V. Oxford University Press. New York.
- Knapp W, *et al.* 1989. Leucocyte Typing IV Oxford University Press. New York.
- McMichael A, *et al.* Eds. 1987. Leucocyte Typing III Oxford University Press. New York.
- Vainer B, *et al.* 2000. *Am. J. Surg. Pathol.* 24:1115. (IHC)
- Ottonello L, *et al.* 1999. *Blood* 93:3505.

6. Metelitsa LS, *et al.* 2002. *Blood* 99:4166.
7. Sadhu C, *et al.* 2007. *J. Leukoc. Biol.* doi:10.1189/jlb.1106680. [PubMed](#)
8. Ihanus E, *et al.* 2007. *Blood* 109:802-810.
9. Gurer C, *et al.* 2008. *Blood* 112:1231. [PubMed](#)
10. Asai A, *et al.* 2009. *J. Lipid Res.* 50:95. [PubMed](#)
11. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
12. Sadhu C, *et al.* 2008. *J. Immunoass. Immunoch.* 29:42. (FC)

Product Citations

1. Buffone A, *et al.* 2018. *J Cell Sci.* 131:. [PubMed](#)
2. Osorio FG *et al.* 2018. *Cell reports.* 25(9):2308-2316 . [PubMed](#)
3. Singh N, *et al.* 2017. *J Endocrinol.* 235:69. [PubMed](#)
4. Matsuyama H, *et al.* 2019. *Sci Rep.* 9:13181. [PubMed](#)
5. Denda-Nagai K, *et al.* 2010. *J Biol Chem.* 285:19193. [PubMed](#)
6. Shen L, *et al.* 2017. *J Virol.* 10.1128/JVI.02310-16. [PubMed](#)
7. Danzer H, *et al.* 2020. *Elife.* 9:00. [PubMed](#)
8. Marjanovic ND, *et al.* 2020. *Cancer Cell.* 38(2):229-246.e13. [PubMed](#)
9. Chakraborty M, *et al.* 2021. *Cell Reports.* 34(2):108609. [PubMed](#)

RRID

AB_314173 (BioLegend Cat. No. 301603)
 AB_314174 (BioLegend Cat. No. 301604)

Antigen Details

Structure	Integrin, type I transmembrane glycoprotein, associates with integrin β_2 (CD18), 145-150 kD
Distribution	Myeloid, dendritic cells, NK cells, B cells and T cell subsets
Function	Adhesion, CTL killing
Ligand/Receptor	CD54, fibrinogen, iC3b, ICAM-1, ICAM-4
Cell Type	Dendritic cells, NK cells, B cells, T cells, Neutrophils, Tregs
Biology Area	Cell Adhesion, Cell Biology, Immunology, Innate Immunity, Neuroscience, Neuroscience Cell Markers, Costimulatory Molecules
Molecular Family	Adhesion Molecules, CD Molecules
Antigen References	<ol style="list-style-type: none"> 1. Petty H. 1996. <i>Immunol. Today</i> 17:209. 2. Springer T. 1994. <i>Cell</i> 76:301. 3. Ihanus E, <i>et al.</i> 2007. <i>Blood</i> 109:802-810.
Gene ID	3687


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

FITC anti-human CD11c, PE anti-human CD11c, Purified anti-human CD11c, PE/Cyanine7 anti-human CD11c, PE/Cyanine5 anti-human CD11c, Biotin anti-human CD11c, APC anti-human CD11c, Alexa Fluor® 488 anti-human CD11c, Alexa Fluor® 647 anti-human CD11c, Pacific Blue™ anti-human CD11c, PerCP/Cyanine5.5 anti-human CD11c, Brilliant Violet 421™ anti-human CD11c, Brilliant Violet 711™ anti-human CD11c, Ultra-LEAF™ Purified anti-human CD11c, Brilliant Violet 510™ anti-human CD11c, Brilliant Violet 605™ anti-human CD11c, Brilliant Violet 650™ anti-human CD11c, Purified anti-human CD11c (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD11c, Brilliant Violet 785™ anti-human CD11c, Alexa Fluor® 700 anti-human CD11c, APC/Fire™ 750 anti-human CD11c

Product Data

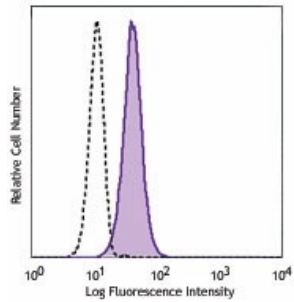

Human peripheral blood granulocytes were stained with CD11c (clone 3.9) FITC (filled histogram) or mouse IgG1, κ FITC isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587